

WONDER PRO FILL

THE STRENGTH & BEAUTY NANO RESTORATIVE

Shade selection in composite restorations is often made confusing by the complex shade recipes advocated by different composite brands. Besides requiring to learn the technique, these also need large inventories of shades & opacities to follow.

Wonder Pro Fill's unique **Bridge Shade Concept** incorporates up to 4 VITA shades in a single syringe; making shade selection easier than ever before. It's best in class filler load (83%) offers exceptional strength and excellent handling.



In addition to making the shade recipe simpler, Wonder Pro Fill incorporates all the required shades and opacities in just 9 syringes. Thus, allowing you to craft aesthetic life-like restorations with minimal inventory.

Best in class **filler load 83%** gives exceptional **compressive strength (330 MPa)** and high abrasion resistance

Exceptional handling - non-sticky, packable consistency makes it a pleasure to work with

Easy and predictable shade match with Bridge shade concept

No shade shift – pre-cure to post cure

Minimal polymerization shrinkage reduces risk of microleakage & cusp fractures

Superior **polishability and polish retention**

With Wonder Pro Fill you can achieve maximum possibilities with minimum inventories.


SIMPLE SHADE RECIPE WITH BRIDGE SHADES

Bridge shade concept with **3 opacities - Body, Dentin & Enamel** enables high end aesthetic restorations with any layering technique - single, dual or multilayering.

VITA* Classic Shade	Wonder Pro Fill Body Shade	Wonder Pro Fill Dentine Shade	Wonder Pro Fill Enamel Shade
B1	A1	0A2	UE
A1			
B2			
D2	A2		
A2			
C2			
D4	A3		
A3			
D3			
B3	A3.5	0A3	
A3.5			
B4			
C3	A4		
A4			
C4			
Layering Technique	Mono		
	Dual Layer		
	Multi Layer		

*VITA is a trademark of VITA Zahnfabrik, Germany

SPECIALLY INDICATED FOR

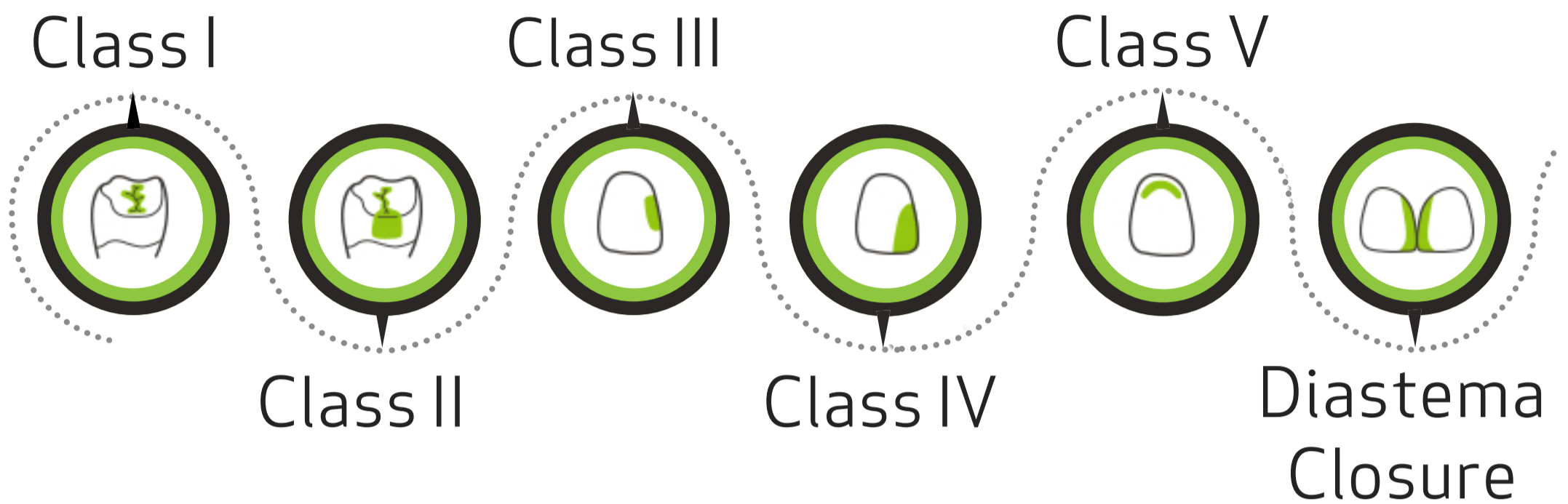

ORDERING INFORMATION

•• 9 Shades:

Body Shades

A1 (A1 → B1)

A2 (B2 → A2 → D2)

A3 (C2 → D4 → D3)

A3.5 (B3 → A3.5 → B4)

A4 (C3 → A4 → C4)

Dentine Shades

0A2 (0A2 → 0B1 → 0B2)

0A3 (0A3 → 0B3 → 0A3.5)

Enamel Shade

UE (UE → CE → NC)

Bleach Shade

XW (XW → W)

 4g Syringe x1

[Click here](#) to find more on what other dentists are talking about Wizdent Products

 fb.com/fineartofdentistry

 www.wizdent.com